



PRODUCT DATA SHEET



XEPTAGEN SpA
VEGA Science Park – Building Auriga
Via delle Industrie, 9
30175 Marghera (VE) – ITALY

Phone: +39 041 509 3910
Fax: +39 041 509 3884

info@xeptagen.com
www.xeptagen.com

© 2004 XEPTAGEN S.p.A.



PRODUCT DATA SHEET

Prostate-IC

XG007

ELISA Kit for the detection of Prostate Antigen (PSA) Immune Complexes (PSA-IgM) in Prostate Cancer

Rev. 04/2010

INTENDED USE

Prostate-IC is an enzyme linked immunosorbent assay (ELISA) for the quantitative measurement of Prostate Specific Antigen (PSA) immune complexes (PSA-IgM) in serum samples.

SUMMARY AND EXPLANATION OF THE TEST

Prostate-IC belongs to a novel generation of in vitro diagnostic methods based on the detection of PSA as circulating Immune Complexes. Prostate-IC is a highly specific and sensitive ELISA for Prostate Cancer detection designed to measure PSA-IgM in patient sera. The amount of PSA-IgM is expressed in Arbitrary Units (AU/mL), using a specific calibrator as reference. Studies have demonstrated that the measurement of PSA-IgM levels in Prostate Cancer patients detects cancer with higher sensitivity with respect to PSA, the well-established serum cancer biomarker, without compromising specificity (1-3). Furthermore, the occurrence of the free and IgM-complexed form of the circulating PSA does not overlap, indicating that PSA-IgM is complementary to PSA and the use of both markers can increase the diagnostic accuracy (1-5).

PRINCIPLE OF THE TEST

Standard Calibrators and specimens are simultaneously incubated with anti-PSA antibodies coated to the wells of a microtiter plate. The immune complexes PSA-IgM are detected by the addition of an enzyme conjugated secondary antibody and an enzyme substrate (TMB). The developed color is proportional to the amount of the analyte in the sample.

REAGENTS AND MATERIALS PROVIDED

XG007-PL: 96 wells multi-strip Assay-Plate, pre-coated with affinity purified rabbit anti-PSA.

XG007-Calibrator: Two vials of calibrator lyophilized from PBS. White powder. Exact concentration on label. Totally soluble.

XG-EA: 200 µL of Enzyme-conjugated goat anti-human IgM secondary antibody (Green cap) 300-fold concentrate in stabilizer solution. The solution contains Proclin as preservative.

XG-CH3: 11 mL of TMB (3,3',5,5'-Tetramethylbenzidine) chromogen solution ready to use.

XG-ST3: 11 mL of 1N HCl Stop solution ready to use.

XG-DB5: 10 mL of concentrated Dilution Buffer solution 5X. Once diluted, the working solution contains 1% BSA

and 0.05% Tween 20 in PBS. The solution contains Proclin as preservative.

XG-WB2: Two tablets of lyophilized Washing Buffer. Once diluted, the working solution contains 0.05% Tween 20 in PBS. Totally soluble.

MATERIAL AND EQUIPMENT REQUIRED

Precision pipettes with disposable tips
Microplate washer
Microplate readers with a 450 ± 20 nm filter
Deionized water

STORAGE CONDITIONS

Storage at 4°C:

XG007-PL, XG-CH3, XG-ST3, XG-DB5*, XG-WB2*

Storage at -20°C:

XG-EA, XG007-Calibrator[§].

Avoid repeated freezing and thawing cycles

(*) Must be used within one month of reconstitution

(§) Must be reconstituted just before the use

EXPIRATION DATE

Expiration date printed on the kit indicates limits of stability.

WARNINGS - POTENTIAL BIOHAZARDOUS MATERIALS

The standard calibrator **XG007-Calibrator** is of **human origin**. The reference material was tested using an approved method of evaluation for the presence of the antibodies to HIV, antibodies to the hepatitis C virus and hepatitis B surface antigens, and found to be negative. **Since no test method can offer complete assurance that HIV, hepatitis B virus, hepatitis C virus, or other infectious agents are absent, all human sourced materials should be considered potentially infectious.** It is recommended that these reagents and human specimens be handled in accordance with the OSHA Standard on Bloodborne Pathogens (6). Biosafety Level 2 (7) or other appropriate biosafety practices (8, 9) should be used for materials that contain or are suspected of containing infectious agents.

SPECIMEN COLLECTION AND PREPARATION

The use of serum samples are recommended for the Prostate-IC assay. Serum specimens should be collected aseptically, avoiding hemolysis if possible.



PRODUCT DATA SHEET

Specimens should be stored at 2-8°C if the assay will be performed within 24 hours after collection. Specimens should be stored frozen if testing will occur after 24 hours.

If frozen, specimens should be mixed thoroughly after thawing to ensure consistency in the results. Avoid repeated freezing and thawing. Specimens showing particulate matter, erythrocytes, or turbidity must be clarified by centrifugation before testing.

INSTRUCTIONS FOR USE

PROCEDURAL NOTES

- Allow samples and reagents to reach room temperature prior to testing. Do not use water baths to thaw samples or reagents.
- Mix samples and all reagents thoroughly before use.
- Avoid excessive foaming of reagents. Also avoid exposure of reagents to excessive heat or light during storage and incubation.
- Avoid handling the tops of the wells both before and after filling.
- Standards and samples should be assayed in duplicate.
- Run a separate standard curve for each assay.
- Use only coated wells from the same reagent batch for each assay. Also do not mix reagents from different kit lots.
- Perform incubations in a sealed box containing a wet paper towel in order to prevent evaporation.

REAGENTS PREPARATION

- Reconstitute XG007-Calibrator with 440 µL of deionized water for each calibrator vial.
- Prepare the required amount of XG-DB5 dilution buffer by diluting 5-fold the concentrated solution in deionized water. If crystals appear upon refrigeration, warm the bottle to 37°C with mixing to dissolve.
- Reconstitute 1 tablet of XG-WB2 washing buffer in 500 mL of deionized water. Prepare 1 liter with two tablets of XG-WB2 washing buffer for the entire plate.
- Prepare the required amount of XG-EA enzyme-conjugated secondary antibody solution diluting 300-fold in reconstituted XG-DB5 dilution buffer.

ASSAY PROTOCOL

1. Prepare assay reagents as described above.
2. Set up the microtiter plate with sufficient wells to enable the running of all required standards and samples.

3. Remove excess microtiter plate strips from the frame and store in the re-sealable foil bag with the desiccant provided.
4. Wash the microtiter plate strips three times with XG-WB2 washing buffer (300 µL/well).
5. Dispense 100 µL/well of standard calibrators (in duplicate) starting from the reconstituted solution and performing in-plate 2-fold serial dilutions in order to obtain a seven-point calibration curve. Use XG-DB5 dilution buffer as diluent (videoclip available at <http://www.xeptagen.com/elisains>). For exact concentration of the reconstituted calibrator please refer to the concentration value (AU/mL) indicated on the XG007-Calibrator vial. Also dispense 100 µL/well of XG-DB5 dilution buffer as blank, in duplicate.
6. Dispense 100 µL/well of a 50-fold diluted sample (in duplicate). Use XG-DB5 dilution buffer as diluent.
7. Incubate 1 hour at room temperature.
8. Wash six times with XG-WB2 washing buffer (300 µL/well).
9. Add 100 µL/well of diluted XG-EA enzyme-conjugated secondary antibody solution.
10. Incubate 1 hour at room temperature.
11. Wash six times with XG-WB2 washing buffer (300 µL/well).
12. Apply 100 µL/well of XG-CH3 chromogen solution.
13. Allow color to develop for 10-15 min at room temperature in the dark then apply 100 µL/well of XG-ST3 Stop Solution and measure OD values of each well using an ELISA plate reader equipped with a 450 nm filter. Stopped reaction should be read within 1 hour.
14. Elaborate OD values with Xerepro software (<http://www.xeptagen.com/software>) or plot the standard curve ΔOD values as described in the next section: Processing of the results.

PROCESSING OF THE RESULTS

Average the duplicate readings for each standard calibrator and sample and subtract the zero standard optical density.

The standard calibrators may be used to construct two distinct standard curves with values reported in AU/mL. Plot on a semi-logarithmic graph the OD readings corresponding to the reconstituted standard calibrator and those corresponding to 1:2 and 1:4 dilutions. Plot on a linear graph the OD readings corresponding to the PSA-IgM titers obtained diluting 1:4 to 1:64 the reconstituted standard calibrator, and the point (0,0).

Data deduction may be performed through computer methods using curve fitting routines or may also be manually deduced on paper.

The immune complexes concentration in the biological sample can be calculated directly from the appropriate



PRODUCT DATA SHEET

standard curve, depending on the sample absorbance value, by interpolation. The value obtained has to be multiplied by the dilution factor of the sample. Samples with OD values exceeding the upper calibration limit should be further diluted and re-measured.

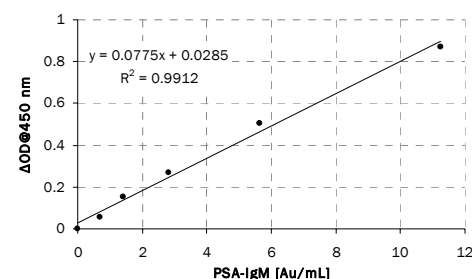


FIG. 1A: range of linearity (0 to 11.25 AU/mL) of a typical standard curve for PSA-IgM after 15 minutes of substrate incubation at room temperature and addition of stop solution

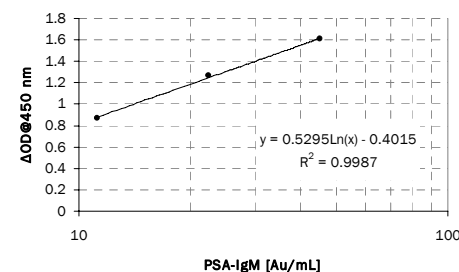


FIG. 1B: range of linearity (from 11.25 to 45.0 AU/mL) of a typical semi-logarithmic standard curve for PSA-IgM after 15 minutes of substrate incubation at room temperature and addition of stop solution

QUALITY CONTROL

The intra- and inter-assay coefficients of variation were determined on 4 typical standard curves and the results were less than 10%.

For optimal performance, the absorbance of the zero standard should be < 0.2 OD450.

It is recommended that each laboratory assays appropriate quality control samples in each run to ensure that all reagents and procedures are correct.

INTERPRETATION

The PSA-IgM cut-off value was 145 AU/mL for differentiating Prostate Cancer from Hyperplasia (2,3).

SPECIFIC PERFORMANCE CHARACTERISTICS

The linear range of the assay is 0-45 AU/mL. The sample with values above 45 AU/mL should be further diluted and re-measured.

REFERENCES

1. Beneduce L, Prayer-Galetti T, Marcello Grimani Giustinian A, Gallotta A, Betto G, Pagano F, Fassina G. **Detection of prostate-specific antigen coupled to immunoglobulin M in prostate cancer patients.** Cancer Detection and Prevention 31: 402-407, 2007.
2. Zani D, Pezzoni G, Pettenò A, Simeone C, Beneduce L, Gallotta A, Fassina G, Cosciani Cunico S. **Validazione dell'uso di PSA e immuno-complessi PSA-IgM nella diagnostica della malattia prostatica organo confinata.** Congresso Nazionale (SIU), Roma, 22-28 September 2008.
3. Prayer-Galetti T, Beneduce L, Dal Moro F, Betto G, Fassina G, Pagano F. **Influenza dell'età dei pazienti, dei valori di PSA sierico e del grado di differenziazione tumorale sui livelli circolanti degli immuno-complessi IgM.** 80° Congresso Nazionale (SIU), Bari, September 27-October 1 2007.
4. Prayer-Galetti T, Beneduce L, Dal Moro F, Betto G, Fassina G, Pagano F. **PSA-IgM complessato (IC) un nuovo biomarcatore per il carcinoma della prostata.** 79° Congresso Nazionale (SIU), Bologna, Italy, June 17-21 2006.
5. Beneduce L, Prayer-Galetti T, Marcello Griman Giustinian A, Gallotta A, Fracalanza S, Betto G, Artibani W, Pagano F, Fassina G. **The occurrence of Prostate Specific Antigen-IgM immune complexes (IC) as novel serum biomarker for prostate cancer.** 21st Annual Congress of European Association of Urology, Paris, 05-08 April 2006. European Urology Supplements 5 (2): 163, 2006.
6. US Department of Labor, Occupational Safety and Health Administration, 29 CFR Part 1910.1030, **Occupational Exposure to Bloodborne Pathogens Final Rule.** Federal Register; 56 (235):64175-82, 1991.
7. US Department of Health and Human Services. **Biosafety in Microbiological and Biomedical Laboratories.** HHS Publication No. (CDC) 93-8395. Washington, DC: US Government Printing Office, May 1993.
8. World Health Organization **Laboratory Biosafety Manual.** Geneva: World Health Organization, 1993.
9. National Committee for Clinical Laboratory Standards. **Protection of Laboratory Workers from Infectious Diseases Transmitted by Blood, Body Fluids and Tissue: Tentative Guideline.** NCCLS Document M29-T2 Villanova: NCCLS: 1-43, 1991.